



www.eyenetra.com

FIND, ALIGN AND CONFIRM

Repeat 3 easy movements to find your refraction results. Grip the NETRA with both hands, wrapping your palms underneath for greater support.


1

Turn the Brightness knob to locate the red and green umbrellas. Make sure both umbrellas are as full as possible before moving on.



2

Turn the Alignment knob to align the two umbrella stems so they overlap completely.

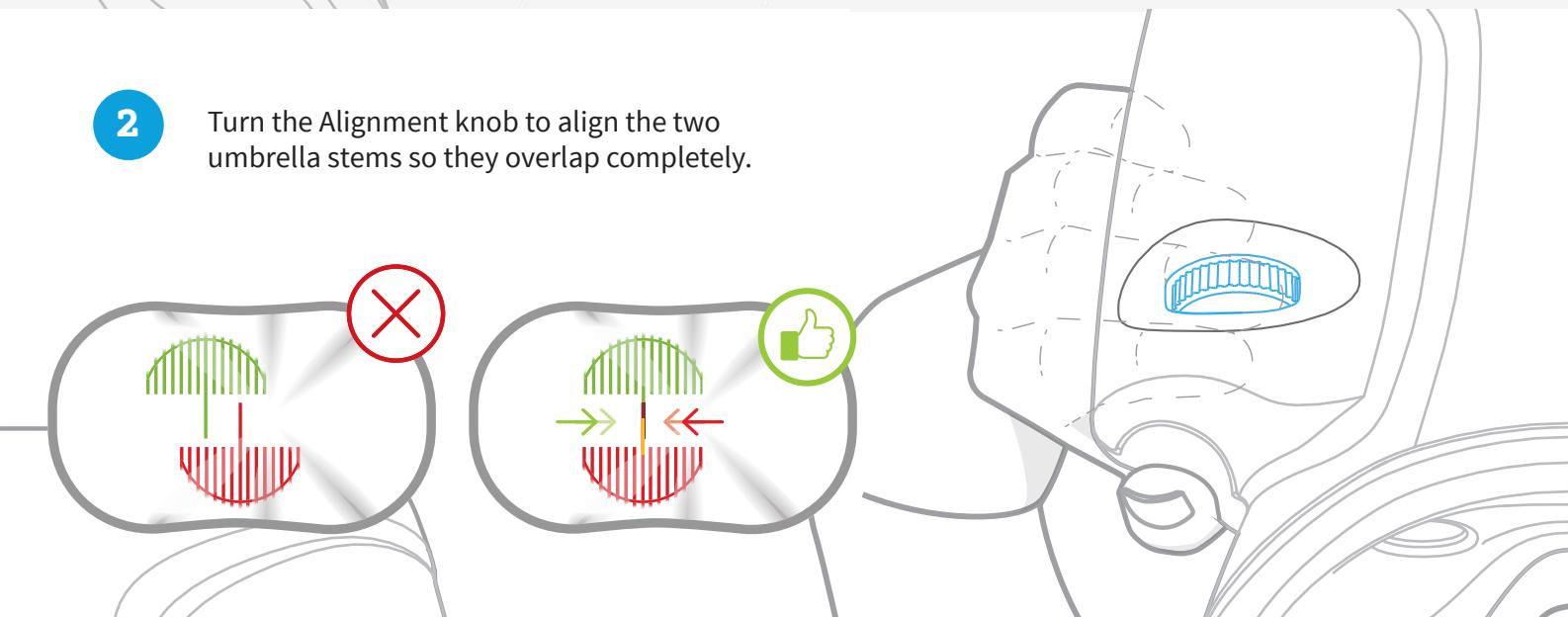

3

Press the Confirm button (✓).



Repeat steps 1, 2 and 3, 8X per eye until you see a blue check mark on screen.

